

#### **UK Office**

#### **Everest Biotech Ltd**

Cherwell Innovation Centre

77 Heyford Park Upper Heyford Oxfordshire OX25 5HD

UK

**Enquiries:** 

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

#### **US Office**

#### **Everest Biotech c/o Abcore**

405 Maple Street, Suite A106 Ramona, CA 92065

USA

Inquiries:

info@everestbiotech.com

Sales:

 $\underline{usasales@everest biotech.com}$ 

Tech support:

support@everestbiotech.com

Tel: 888-320-4628 (toll-free) Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

### EB06377-T - Goat Anti-GAPDH (C Terminus) Loading Control Antibody - Trial

Size: 20µg specific antibody in 40µl



### **Target Protein**

**Principal Names:** GAPDH, glyceraldehyde-3-phosphate dehydrogenase, HGNC:4141, G3PD, GAPD, MGC88685, aging-associated gene 9 protein, glyceraldehyde 3-phosphate dehydrogenase, epididymis secretory sperm binding protein Li 162eP, HEL-S-162eP,

peptidyl-cysteine S-nitrosylase GAPDH

Official Symbol: GAPDH

Accession Number(s): NP\_002037.2

Human GenelD(s): 2597

**Important Comments:** GAPDH is constitutively expressed in almost all tissues at high levels. It is therefore a useful marker when a loading/positive control is required in western blotting.

#### **Immunogen**

Peptide with sequence C-HQVVSSDFNSDT, from the C Terminus of the protein sequence according to NP\_002037.2.

Please note the peptide is available for sale.

### **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

#### **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:1000.

**Western blot:** Approx. 37kDa band observed in lysates of cell lines HeLa, HEK293, LNCaP and A549, and in Human Liver lysates (calculated MW of 36.1kDa according to NP\_002037). Recommended concentration: 0.1-0.5μg/ml. Primary incubation 1 hour at room temperature.

**Positive Control:** A batch specific positive control lysate is available for this product. Please contact Sales@everestbiotech.com for availability.

**IHC:** Paraffin embedded Human Kidney and Lung. Recommended concentration: 5-7µg/ml.

**Immunofluorescence:** Strong expression of the protein seen in the cytoplasm and plasma membranes of HeLa cells. Recommended concentration: 10µg/ml.

#### **Species Reactivity**

Tested: Human

Expected from sequence similarity: Human, Mouse, Dog, Pig

#### **Specific References**

This antibody (previous batch) has been successfully used in Western blot on Drosophila:

George J, Jacobs HT

Germline knockdown of spargel (PGC-1) produces embryonic lethality in Drosophila.

Mitochondrion. 2019 Aug 29;49:189-199

PMID: 31473309

# This antibody (previous batch) has been successfully used in Western blot on Drosophila:

Jack George & Howard T. Jacobs

Minimal effects of spargel (PGC- $1\alpha$ ) overexpression in a Drosophila mitochondrial disease model

bioRxiv 529545; doi: https://doi.org/10.1101/529545 (2019)

PMID: 31292108

# This antibody (previous batch) has been successfully used in Western blot on Human:

Ananda Ayyappan Jaguva Vasudevan, Ulrike Kreimer, Wolfgang A Schulz, Aikaterini Krikoni, Gerald G Schumann, Dieter Häussinger, Carsten Münk, Wolfgang Goering APOBEC3B Activity Is Prevalent in Urothelial Carcinoma Cells and Only Slightly Affected by LINE-1 Expression

Front Microbiol. 2018 Sep 4;9:2088.

PMID: 30233553

# This antibody (previous batch) has been successfully used in Western blot on Drosophila:

Cagri Yalgin, Bohdana Rovenko, Ana Andjelkovi**■**, Margot Neefjes, Burak Oymak, Eric Dufour, Ville Hietakangas & Howard T. Jacobs

Effects on dopaminergic neurons are secondary in COX-deficient locomotor dysfunction in Drosophila

https://ssrn.com/abstract=3600553 (3rd June 2020)

PMID: 0

#### This antibody (previous batch) has been successfully used in the following paper:

Krzysztof Sikorski, Adi Mehta, Marit Inngjerdingen, Flourina Thakor, Simon Kling, Tomas Kalina, Tuula A. Nyman, Maria Ekman Stensland, Wei Zhou, Gustavo A. De Souza, Lars Holden, Jan Stuchly, Markus Templin and Fridtjof Lund-Johansen

A high-throughput pipeline for validation of antibodies

Nat Methods. 2018 Nov;15(11):909-912

PMID: 30377371

# This antibody (previous batch) has been successfully used in Western blot on Human:

Osei Kuffour E, Schott K, Jaguva Vasudevan AA, Holler J, Schulz WA, Lang PA, Lang KS, Kim B, Häussinger D, König R, Münk C

USP18 (UBP43) Abrogates p21-Mediated Inhibition of HIV-1

J Virol. 2018 Sep 26;92(20). pii: e00592-18

PMID: 30068654

# This antibody (previous batch) has been successfully used in Western blot on Human and Mouse:

Jaguva Vasudevan AA, Bähr A, Grothmann R, Singer A, Häussinger D, Zimmermann A, Münk C.

MXB inhibits murine cytomegalovirus

Virology. 2018 Sep;522:158-167.

PMID: 30032029

### This antibody (previous batch) has been successfully used in Western blot on Human:

Hain A, Krämer M, Linka RM, Nakhaei-Rad S, Ahmadian MR, Häussinger D, Borkhardt A, Münk C.

IL-2 Inducible Kinase ITK is Critical for HIV 1 Infection of Jurkat T-cells.

Sci Rep. 2018 Feb 16;8(1):3217.

PMID: 29453458

# This antibody (previous batch) has been successfully used in Western blot on Drosophila:

Scialo F, Sriram A, Stefanatos R, Sanz A.

Practical Recommendations for the Use of

the GeneSwitch Gal4 System to Knock-Down Genes in Drosophila melanogaster.

PLoS One. 2016 Aug 29;11(8)

PMID: 27570965

# This antibody (previous batch) has been successfully used in Western blot on Human:

Happold C, Roth P, Silginer M, Florea AM, Lamszus K, Frei K, Deenen R, Reifenberger G, Weller M.

Interferon- $\beta$  induces loss of spherogenicity and overcomes therapy resistance of glioblastoma stem cells.

Mol Cancer Ther. 2014 Apr;13(4):948-61.

PMID: 24526161

# This antibody (previous batch) has been successfully used in Western blot on Human:

Happold C, Roth P, Wick W, Schmidt N, Florea AM, Silginer M, Reifenberger G, Weller M. Distinct molecular mechanisms of acquired resistance to temozolomide in glioblastoma cells

J Neurochem. 2012 May 7. doi: 10.1111/j.1471-4159.2012.07781.x.


PMID: 22564186

### This antibody (previous batch) has been successfully used in Western blot on Rat:


Kiepe D, Van Der Pas A, Ciarmatori S, Ständker L, Schütt B, Hoeflich A, Hügel U, Oh J, Tönshoff B.

Defined carboxy-terminal fragments of insulin-like growth factor (IGF) binding protein 2 exert similar mitogenic activity on cultured rat growth plate chondrocytes as IGF-I. Endocrinology. 2008 Oct;149(10):4901-11.


PMID: 18556354

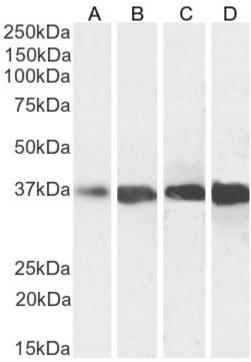


EB06377 (0.5μg/ml) staining of HeLa (A), HEK293 (B), LNCaP (C) and A549 (D) cell lysate (35μg protein in RIPA buffer). Detected by chemiluminescence.

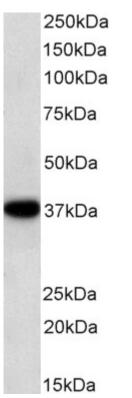


EB06377 (0.5μg/ml) staining of Human Liver lysate (35μg protein in RIPA buffer). Detected by chemiluminescence.

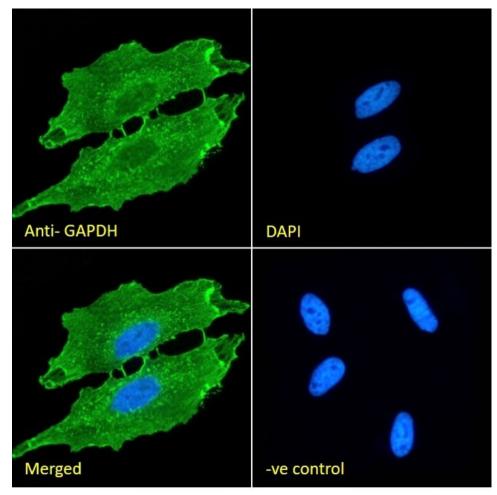


EB06377 Immunofluorescence analysis of paraformaldehyde fixed HeLa cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing membrane and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

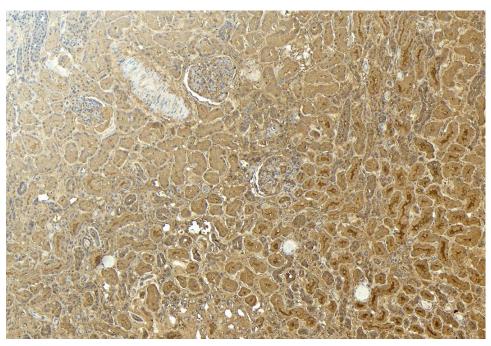


EB06377 (5µg/ml) staining of paraffin embedded Human Kidney. Heat induced antigen retrieval with citrate buffer pH 6, HRP-staining.

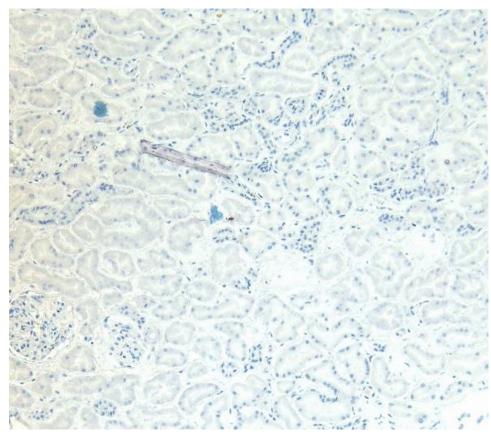

 ${\tt EB06377\ Negative\ Control\ showing\ staining\ of\ paraffin\ embedded\ Human\ Kidney,\ with\ no\ primary\ antibody.}$ 

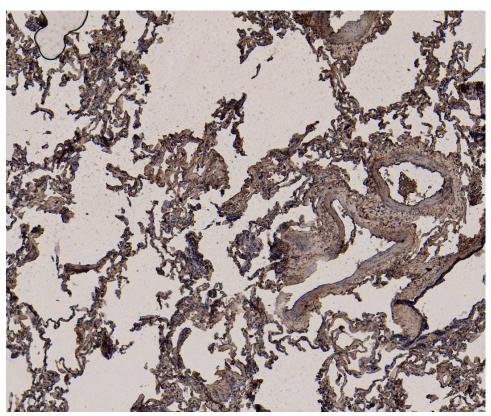

EB06377 ( $6\mu g/ml$ ) staining of paraffin embedded Human Lung. Heat induced antigen retrieval with citrate buffer pH 6, HRP-staining.

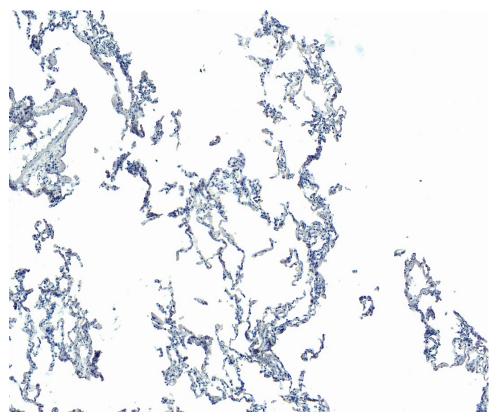

 ${\sf EB06377\ Negative\ Control\ showing\ staining\ of\ paraffin\ embedded\ Human\ Lung,\ with\ no\ primary\ antibody.}$